Botón Septal

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PJ SMITH MEDICAL PRODUCTS MARKETING
18 Yeates Close • Thame OX93AR, UK

INVOTEC INTERNATIONAL, INC.

6833 Phillips Industrial Blvd. • Jacksonville, Florida 32256 U.S.A.
Tel: 800-998-8580 / 904-880-1229 • Fax: 904-886-9517 • www.invotec.net

PJ SMITH MEDICAL PRODUCTS MARKETING
18 Yeates Close • Thame OX93AR, UK

Septal Button

Please Read Carefully:

NOTE: THE MEDICAL PROFESSIONAL IS RESPONSIBLE FOR PROPER SURGICAL PROCEDURES AND TECHNIQUES. THE FOLLOWING INSTRUCTIONS ARE FOR INFORMATION ONLY. THE ATTENDING SURGEON MUST EVALUATE THE APPROPRIATENESS OF THE PROCEDURE BASED ON HIS/HER OWN MEDICAL TRAINING AND EXPERIENCE.

Indications

When a patient is a poor candidate for surgical correction of a septal perforation, and the nature of the perforation, that is size and location, and septal anatomy are appropriate, closure with a septal button is indicated. Generally, anterior septal perforations that are 2cm. or less in maximum diameter are suitable for closure with a septal button.

Contraindications

Septal deformities are relative contraindications for septal button placement, as the flange on the convex side of the deformity may increase nasal obstruction.

How to use the Septal Button

Step 1: After adequate cocainization of the nose, 1% lidocaine with 1:100,000 epinephrine is infiltrated about the margins of the perforation.

Step 2: Trim the flanges of the button to permit 3-4mm. of overlap with respect to the margins of the perforation.

Step 3: Choose one of the following two methods of insertion:

METHOD 1:

Collapse one flange and introduce the button into the perforation, then gently pull it through from the opposite end.

METHOD 2:

Sew the edges of one flange together in a rosette manner with a suture, tie suture and leave long. Pass the long end of the suture through the perforation and gently guide the button in place using a hemostat from the back side to apply gentle pressure. Cut and remove the suture allowing the sewn flange to open flat.

NOTE: CARE MUST BE TAKEN SO THAT THE FLANGES OF THE BUTTON DO NOT ABUTT THE JUNCTION OF THE SEPTUM WITH THE UPPER LATERAL CARTILAGE, AS THIS IS A PARTICULARLY ANNOYING SENSATION FOR THE PATIENT.

Removal

Should it become necessary to remove the button, it is best removed by carefully cutting the center post with a scalpel and removing each piece of the button through separate nostrils.

Caution: U.S. Federal law restricts this device to sale by or on the order of a physician.

Warranty


Invotec International, Inc. warrants that the product is free from defects in material and workmanship. Invotec will replace or provide a refund for any product found to be defective so long as the product is returned according to the Returned Goods instructions in the Sales Policy. Invotec shall not be liable for any consequential loss, damage or expense directly or indirectly arising from the use of, or inability to use, this product. THE FOREGOING WARRANTY IS IN LIEU OF AND EXCLUDES ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, HOWEVER ARISING, INCLUDING MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, AGAINST INFRINGEMENT OR OTHERWISE.

Invotec International neither assumes, nor authorizes any person to assume for it, any other additional liability or responsibility with respect to this product.

INVOTEC INTERNATIONAL, INC.

6833 Phillips Industrial Blvd. • Jacksonville, Florida 32256 U.S.A.
Tel: 800-998-8580 / 904-880-1229 • Fax: 904-886-9517 • www.invotec.net

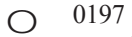PJ SMITH MEDICAL PRODUCTS MARKETING
18 Yeates Close • Thame OX93AR, UK

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